

GE Healthcare

illustra

GenomiPhi HY DNA Amplification Kit

Product Web Protocol

Codes: 25-6600-20
25-6600-22
25-6600-25



Contents

1. Legal	3
1.1. Product use restriction/warranty	3
2. Handling, preparation, and storage of starting materials	4
2.2. Storage Conditions	4
2.3. Expiry	4
3. Components	5
3.1. Kit Contents	5
3.2. Reagents to be supplied by the user	5
4. Product Description	6
4.1. The Basic Principle	6
4.2. Kit Specifications	6
5. Protocols	8
5.1. Short Protocols	8
5.2. General Protocol	8
5.3. Chemical denaturation of template DNA	9
5.4. Quantitation of Amplification Products	11
6. Appendix	14
6.1. Troubleshooting	14
6.2. References	15
6.3. Related Products	15

1. Legal

GE, imagination at work and GE Monogram are trademarks of General Electric Company.

illustra, Autoseq, CyScribe, DYEnamic, GenomiPhi, GFX, MegaBACE, MicroSpin, ProbeQuant, PureTaq and Ready-To-Go are trademarks of GE Healthcare companies.

Nucleon is a trademark of Tepnel Life Sciences Plc.

PicoGreen is a trademark of Invitrogen, Inc.

For use only as licensed by Qiagen GmbH and GE Healthcare Bio-sciences Corp. The Phi 29 DNA polymerase enzyme may not be re-sold or used except in conjunction with the other components of the **GenomiPhi HY DNA Amplification Kit**. See U.S. patent Nos. 5,854,033, 6,124,120, 6,143,495, 5,001,050, 5,198,543, 5,576,204 and related U.S. and foreign patents.

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplied them. A copy of these terms and conditions is available on request.

Contact your GE Representative for the most current information and a copy of the terms and conditions.

© 2006-2007 General Electric Company - All rights reserved.

First published 2006

<http://www.gehealthcare.com/lifesciences>
GE Healthcare UK Limited Amersham Place,
Little Chalfont, Buckinghamshire,
HP7 9NA, UK.

1.1. Product use restriction/warranty

The **GenomiPhi HY DNA Amplification Kit** components have been designed, developed, and sold **for research purposes only**. They are suitable **for *in vitro* use only**. Not recommended or intended for diagnosis of disease in humans or animals. No claim or representation is intended for its use to identify any specific organism or for clinical use (diagnostic, prognostic, therapeutic, or blood banking). Do not use internally or externally in humans or animals.


It is the responsibility of the user to verify the use of the **GenomiPhi HY DNA Amplification Kit** for a specific application range as the performance characteristic of this kit has not been verified to a specific organism.

2. Handling. Safety Warnings and Precautions

2.1. Handling, preparation, and storage of starting materials

 **This kit is sensitive to small amounts of DNA. Wear gloves and safety glasses to avoid contamination.**

This product and its components should be handled only by persons trained in laboratory techniques and used in accordance with the principles of good laboratory practice. All chemicals should be considered potentially hazardous; therefore, when handling chemical reagents, it is advisable that suitable protective clothing, such as laboratory overalls, safety glasses and gloves be worn. Care should be taken to avoid contact with skin or eyes. In case of contact with skin or eyes, wash immediately with water. See the appropriate Material Safety Data Sheet for specific recommendations.

 **GenomiPhi™ HY Amplification Kit is optimized for whole genome amplification from at least 10 ng of high quality genomic DNA template.** Use of less DNA or low quality DNA (such as degraded DNA, or DNA from formalin fixed paraffin embedded samples) can result in amplification bias.

In the absence of template DNA an amplification product will still be produced.

This product is the result of amplification of the DNA hexamers and will not interfere with downstream applications. The yield from no DNA reactions will be the same as from with DNA reactions.

 **Wear gloves at all times during the preparation to avoid contamination.**

2.2. Storage conditions

Store the kit at -70°C.

The enzyme mix must be stored at -70°C; all other components may be stored at -20°C. Thaw components on ice and maintain at 0°C to 4°C during handling.

2.3. Expiry

This product has been designed to deliver high quality results for up to 18 months from the manufacturing date. Please refer to the detail expiration date on the product label.

3. Components

3.1 Kit Contents

GenomiPhi HY DNA Amplification Kit

	25 Reactions	100 Reactions	1000 Reactions
Cat. No.	25-6600-22	25-6600-20	25-6600-25
Sample Buffer (Green Cap)	565 µl	2.25 ml	10 x 2.25 ml
Reaction Buffer (Blue Cap)	565 µl	2.25 ml	10 x 2.25 ml
Enzyme Mix (Yellow Cap)	65 µl	250 µl	10 x 250 µl
Control DNA (Lambda), 10 ng/µl	20 µl	20 µl	10 x 20 µl

3.2. Reagents to be supplied by the user

Liquid-handling supplies - Vials, pipettes, microcentrifuge, and vacuum centrifuge. Perform all amplification reactions in plastic microcentrifuge tubes (typically 0.5 ml), or in 96-well or 384-well plates suitable for sealing and incubating at 30°C.

Thermocycler or water baths - for incubations at 30°C, 65°C and 95°C.

4. Product description

4.1 The basic principle

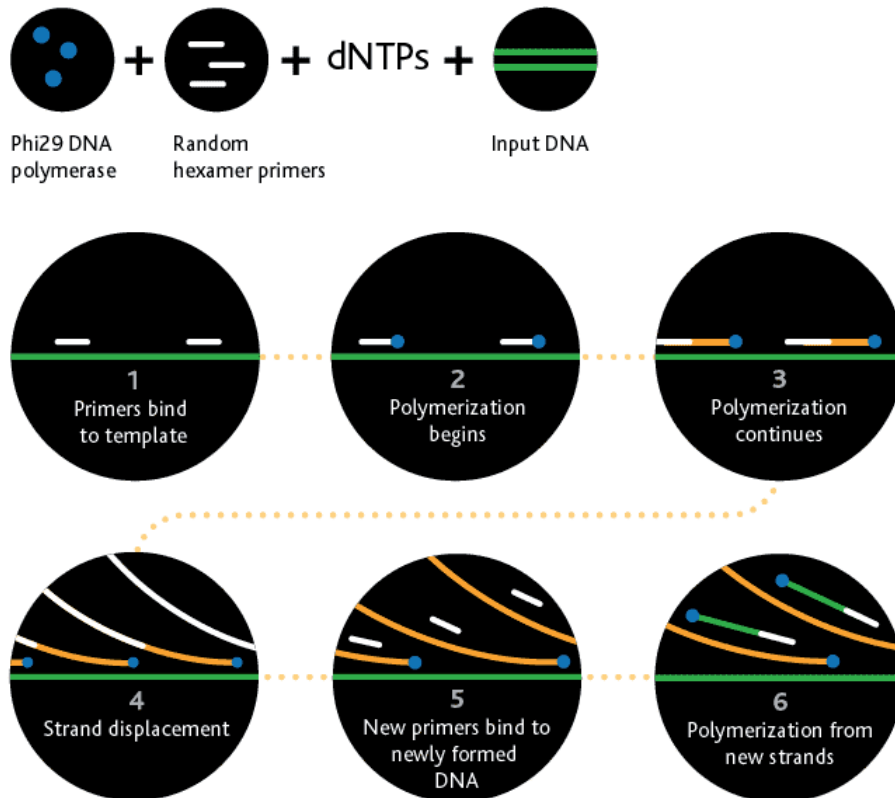


Figure 1 Overview of the **GenomiPhi HY DNA Amplification Kit** procedure

Figure 1 shows an overview of whole genome amplification by isothermal stand displacement using the **GenomiPhi HY DNA Amplification Kit**. DNA is briefly heat-denatured then cooled in sample buffer containing random hexamers that non-specifically bind to the DNA. A master-mix containing DNA polymerase, additional random hexamers, nucleotides, salts and buffers is added and isothermal amplification proceeds at 30°C for 4 hours. After amplification the enzyme is heat inactivated during 10 minute incubation at 65°C.

4.2 Kit specifications

Typical amplification kinetics with **GenomiPhi HY DNA Amplification Kit** is shown in Figure 2 below. Microgram quantities of DNA are generated from nanogram amounts of starting material in only a few hours. Typical DNA yields from a **GenomiPhi HY DNA Amplification Kit** reaction are 45 µg per 50 µl reaction with an average product length of greater than 10 kb. DNA replication is extremely accurate due to the proofreading 3'-5' exonuclease activity of the DNA polymerase (3, 4). Amplification reactions are generally complete in 3-4 h, yielding approximately 45 µg DNA. Reactions containing no DNA will produce the same amount as with DNA reactions.

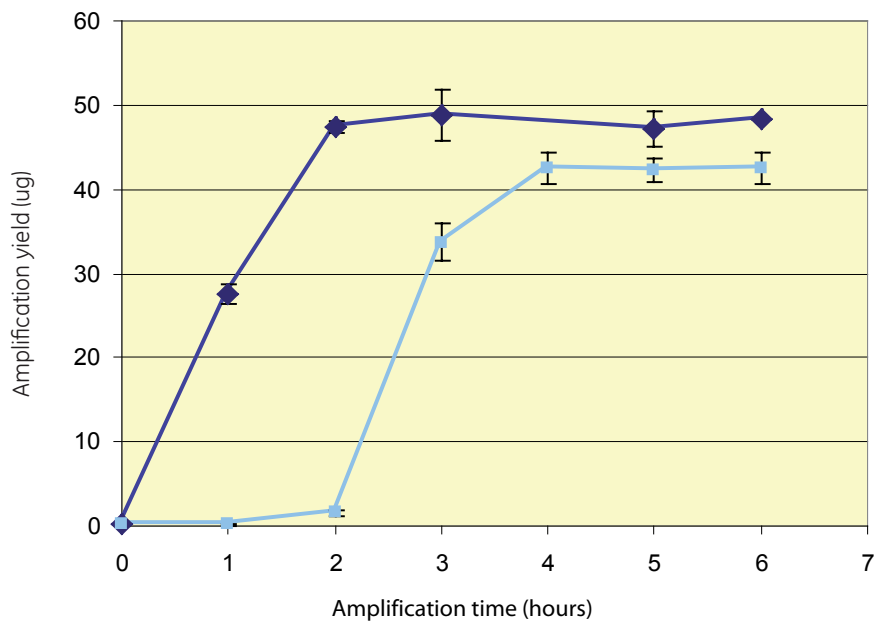


Figure 2 Amplification kinetics of a **GenomiPhi HY DNA Amplification Kit** reaction.

The starting material for **GenomiPhi HY DNA Amplification Kit** reactions can be purified DNA or non-purified cell lysates. Most commercial DNA isolation kits and homemade purification procedures produce suitable DNA for the amplification. Protocols for the amplification of DNA from various clinical samples including blood and buccal cells are available from our web site (www.genomiphi.com).

5. Protocols

5.1. Short Protocol

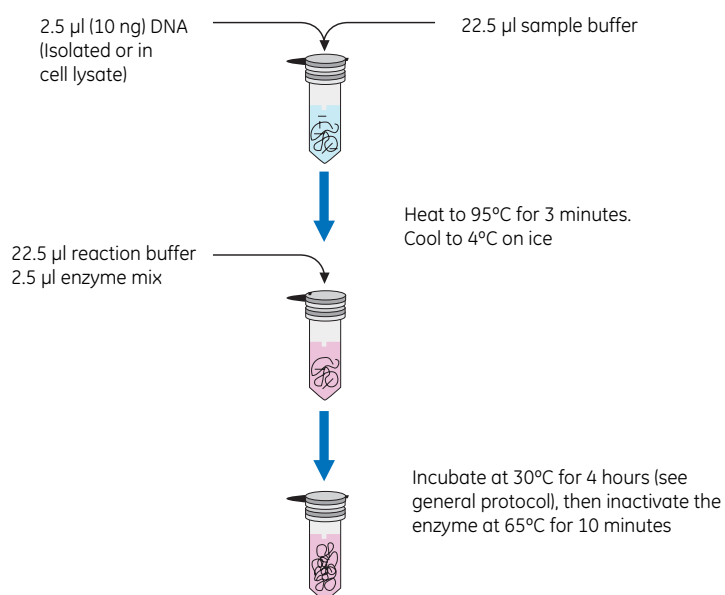


Figure 4 Schematic representation of GenomiPhi HY DNA Amplification Kit protocol.

5.2 General Protocol

The steps outlined below describe a general protocol for amplifying template DNA. This protocol should be considered a starting point for optimizing the reaction in your laboratory.

1 Mix sample buffer with template DNA

Add **22.5 µl** Sample Buffer to **2.5 µl** of 10 ng template DNA.

22.5 µl Sample Buffer

2.5 µl DNA Template

(10 ng)



Template DNA should be resuspended in TE or water.

2 Denature template DNA

Heat the samples to **95°C** for **3 minutes** then cool to **4°C** on ice.

95°C

for 3 minutes



Heating the DNA for longer than 3 minutes or at higher temperatures can cause damage to the DNA.

3 Prepare the master mix for each amplification reaction

For each amplification reaction, on ice combine 22.5 µl of Reaction Buffer with 2.5 µl of Enzyme Mix.



Prepare the master mix only in sufficient quantities and immediately prior to use. Keep the master mix on ice and discard any unused portion. The master mix contains all the components required for DNA amplification and will generate amplification products if exposed to temperatures > 4°C for sufficient time.

4 Transfer master mix to cooled sample.

Add 25 µl of prepared master mix from **Step 3** to each cooled sample, on ice.

25 µl
Master Mix

5 Incubate for DNA amplification.

Incubate the samples at 30°C for 4 hours.

30°C
for 4 hours

6 Phi29 DNA polymerase enzyme.

Heat the samples to 65°C for 10 minutes then cool to 4°C.



Heating is required to inactivate the exonuclease activity of the DNA polymerase which would otherwise begin to degrade the amplification product.

65°C
for 10 minutes

7 Storage of amplified material.

Store amplification reactions at -20°C.



GenomiPhi DNA HY amplification products should be stored and treated as genomic DNA.

Due to the viscosity of the amplification product, dilution in 2 volumes of TE can ease handling. Mix amplification products thoroughly after freezing and thawing.

5.3 Chemical denaturation of template DNA

The steps outlined below describe a general protocol for amplifying template DNA.






1 Prepare denaturation and neutralization buffers.

Denaturation Buffer

400 mM KOH
10 mM EDTA

Neutralization Buffer

200 mM HCl
300 mM Tris-HCl, pH 7.5

2	Denature template DNA Mix 2.5 µl of template (10 ng total) with 2.5 µl of Denaturation Buffer at room temperature . Incubate for 3 minutes.	2.5 µl of DNA + 2.5 µl Denaturation Buffer
	 Mix by pipetting up and down. Do not vortex.	RT for 3 min
3	Neutralize DNA solution. Add 5 µl of Neutralization Buffer.	5 µl Neutralization Buffer
	 Mix by pipetting up and down. Do not vortex.	
4	Prepare the master mix for each amplification reaction For each amplification reaction, on ice combine 15 µl of Sample Buffer and 22.5 µl of Reaction Buffer, with 2.5 µl of Enzyme Mix.	
	 Prepare the master mix only in sufficient quantities and immediately prior to use. Keep the master mix on ice and discard any unused portion. The master mix contains all the components required for DNA amplification and will generate amplification products if exposed to temperatures > 4°C for sufficient time.	
5	Transfer master mix to cooled sample. Add 40 µl of prepared master mix from Step 4 to each cooled sample, on ice.	40 µl Master Mix
6	Incubate for DNA amplification. Incubate the samples at 30°C for 4 hours .	30°C for 4 hours
7	Phi29 DNA polymerase enzyme. Heat the samples to 65°C for 10 minutes then cool to 4°C .	
	 Heating is required to inactivate the exonuclease activity of the DNA polymerase which would otherwise begin to degrade the amplification product.	65°C for 10 minutes
8	Storage of amplified material. Store amplification reactions at -20°C .	
	 GenomiPhi DNA HY amplification products should be stored and treated as genomic DNA.	

Due to the viscosity of the amplification product, dilution in 2 volumes of TE can ease handling. Mix amplification products thoroughly after freezing and thawing.

5.4 Quantitation of amplification products

Quantification is generally not required as every reaction will yield approximately the same amount of DNA. **Quant-iT™ PicoGreen® dsDNA quantitation reagent** (Invitrogen, P7581) is recommended if accurate quantitation is required.



Quantitation of non-purified amplification products by UV absorption will generate inaccurate results due to the presence of unused hexamers in the completed reaction.

1 Prepare TE buffer.

Dilute the concentrated TE buffer included in the kit to 1x concentration using water.

Dilute 1:20



Use only sterile, distilled, DNase free water when preparing the dilution to ensure accurate quantitation.

2 Prepare 1:25 dilution of PicoGreen reagent.

Determine the required volume of a 1:25 dilution of PicoGreen reagent.

- Volume = 100 µl/sample x # of samples

Determine the volume of stock PicoGreen reagent necessary to produce the required volume of a 1:25 dilution

- Volume = $\frac{\text{volume of required dilution}}{25}$



Reagent adsorbs to glass surfaces. Use plasticware only. Protect the solution from light at all times.

3 Prepare the λ DNA standard curve.

Dilute the λ DNA standard supplied in the kit to a **10 ng/µl** working solution. Use this working stock to prepare a standard curve (see below). **Add 100 µl** of each dilution to each well of the assay plate.

Standard Number	λ DNA (10 ng)	λ DNA (10 ng/µl)	1 x TE
1	600	60 µl	40 µl
2	500	50 µl	50 µl
3	400	40 µl	60 µl
4	200	20 µl	80 µl
5	100	10 µl	90 µl
6	50	5 µl	95 µl
7	25	2.5 µl	97.5 µl
8	0	0 µl	100 µl

4 Dilute the GenomiPhi HY amplification products twice.

Dilute the GenomiPhi HY amplification products 1:4 by **adding 150 µl** of 1x TE to each amplification reaction.

Dilute the 1:4 dilution in a separate tube by adding **10 µl** of 1:4 dilution to **90 µl** of 1 x TE.

Serially Dilute
1:4
+
1:10



Due to the viscosity of the amplification product, mix amplification products thoroughly by vortexing heavily.

5 Add diluted GenomiPhi amplification products to the assay plate.

Aliquot 90 µl of 1x TE into each sample well. Add 10 µl of diluted sample for a final volume of 100 µl.

90 µl of TE
+
10 µl of diluted sample



Because the amplification product is diluted **before** the assay, the dilution factor must be taken into consideration when calculating total yields.

6 Add diluted PicoGreen to sample wells.

Add 100 µl of the 1:25 dilution of PicoGreen to all wells containing standards and samples. Mix contents well by pipetting up and down.

Seal the plate with foil and spin in microplate centrifuge for 1 minute at < 200 x g to eliminate bubbles.



Protect plate from light at all times. The plate must be read 5-10 min after addition of PicoGreen reagent to ensure accurate quantitation.

7 Measure the sample fluorescence.

Place the sample assay plate into a fluorescence microplate reader.

Set the fluorescence reader at the following parameters:

- Excitation wavelength: **480 nm**
- Emission wavelength: **520 nm**
- Gain: **Optimal**



If it is not possible to set the instrument gain to optimal, find a way to have the instrument read the sample where the highest DNA concentration generates readings that fall within the linear dynamic range of the instrument.

8 Calculate the concentration of the amplification product.

Correct the data by subtracting the value obtained for the blank from that of each of the samples. Generate a standard curve of fluorescence versus DNA concentration using the corrected data. Determine the concentration of GenomiPhi HY Kit amplified products from the equation of the line derived from the standard curve.

6. Appendix

6.1 Troubleshooting

Problem	Possible cause and suggestions
Reduced yield/ no amplification product	Contamination of template DNA <ul style="list-style-type: none">• Excessive contaminants carried over from the starting material can inhibit the DNA polymerase. Dilute or clean-up the DNA and re-amplify.• Extending the amplification time will help when inhibitory material is causing reduced yields.
	Inactive Enzyme <ul style="list-style-type: none">• It is critical that the enzyme be stored properly. The Enzyme Mix should be stored at -70°C. If the material will be consumed within 2 months, -20°C storage may be used. The freezer must not be a frost-free unit.• Perform a control reaction to confirm performance of the enzyme.
	Low quality DNA <ul style="list-style-type: none">• Amplification kinetics strongly favors intact templates. Avoid template preparation steps that can damage DNA.
	Prolonged denaturation <ul style="list-style-type: none">• Heating at 95°C for 3 minutes is sufficient to denature template DNA and facilitate primer annealing. Longer denaturing times can nick the template and decrease the amplification efficiency.
Poor performance in downstream applications	Degraded/low amounts of template DNA <ul style="list-style-type: none">• In the absence of input DNA, there will most likely still be an amplification product. This amplification product will be nonspecific and may not generate expected results in downstream applications.• Degraded or low amounts of starting DNA template may not amplify consistently or representatively. Increase the amount of starting DNA.• Use high quality genomic DNA for amplification.
	Inhibition of optimized downstream conditions <ul style="list-style-type: none">• Starting material components can inhibit amplification reaction. Purify the starting material using a suitable column prior to amplification.• For some downstream applications, components of the GenomiPhi HY reaction will alter previously optimized downstream conditions. Purify the amplification products using a suitable column after amplification.

6.2 References

1. Dean, F. *et al.*, *Genome Research* **11**, 1095–1099 (2001).
2. Lizardi, P. *et al.*, *Nat. Genet.* **19**, 225–232 (1998).
3. Estaban, J.A. *et al.*, *J. Biol. Chem.* **268**, 2719–2726 (1993).
4. Nelson, J.R. *et al.*; *BioTechniques* **32**, S44–S47 (2002).

6.3 Related Products

DNA Purification Products

GFX™ Genomic Blood DNA Purification Kit
Nucleon™ BACC Genomic DNA Extraction Kits
GFX™ Micro Plasmid Prep Kit
MicroSpin™ Columns
GFX™ PCR DNA and Gel Band Purification Kit
CyScribe™ GFX™ Purification Kit
ProbeQuant™ G-50 Micro Columns

PCR Products

PureTaq™ Ready-To-Go™ PCR Beads
FideliTaq™ PCR Master Mix
dNTP Set, 100 mM Solution (dATP, dGTP, dCTP, dTTP)
Exo-SAP IT
Taq Polymerase

Sequencing Products

MegaBACE™ DNA Analysis System
DYEnamic™ ET Dye Terminator Cycle Sequencing Kit
Autoseq™ G-50

GE Healthcare offices:

GE Healthcare Bio-Sciences AB
Björkgatan 30 751 84
Uppsala Sweden

GE Healthcare Europe GmbH
Munzinger Strasse 5 D-79111
Freiburg Germany

GE Healthcare UK Limited
Amersham Place
Little Chalfont Buckinghamshire
HP7 9NA UK

GE Healthcare Bio-Sciences Corp
800 Centennial Avenue
P.O. Box 1327
Piscataway
NJ 08855-1327 USA

GE Healthcare Bio-Sciences KK
Sanken Bldg. 3-25-1
Hyakunincho Shinjuku-ku
Tokyo 169-0073
Japan

GE Healthcare regional office contact numbers:**Asia Pacific**

Tel: +65 6275 1830
Fax: +65 6275 1829

Australasia

Tel: + 61 2 9899 0999
Fax: +61 2 9899 7511

Austria

Tel: 01/57606-1619
Fax: 01/57606-1627

Belgium

Tel: 0800 73 888
Fax: 02 416 82 06

Canada

Tel: 1 800 463 5800
Fax: 1 800 567 1008

Central, East, & South East Europe

Tel: +43 1 982 3826
Fax: +43 1 985 8327

Denmark

Tel: 45 16 2400
Fax: 45 16 2424

Eire

Tel: 01494 544000
Fax: 0044 1494 542010

Finland & Baltics

Tel: +358-(0)9-512 39 40
Fax: +358 (0)9 512 39 439

France

Tel: 01 6935 6700
Fax: 01 6941 9677

Germany

Tel: (089) 96281 660
Fax: (089) 96281 620

Greater China

Tel: +852 2100 6300
Fax: +852 2100 6338

Italy

Tel: 02 27322 1
Fax: 02 27302 212

Japan

Tel: +81 3 5331 9336
Fax: +81 3 5331 9370

Latin America

Tel: +55 11 3933 7300
Fax: + 55 11 3933 7304

Middle East & Africa

Tel: +30 210 9600 687
Fax: +30 210 9600 693

Netherlands

Tel: 0800 8282821
Fax: 0800 8282824

Norway

Tel: 815 65 555
Fax: 815 65 666

Portugal

Tel: 21 417 7035
Fax: 21 417 3184

Russia & other C.I.S. & N.I.S

Tel: +7 495 956 5177
Fax: +7 495 956 5176

Spain

Tel: 93 594 49 50
Fax: 93 594 49 55

Sweden

Tel: 018 612 1900
Fax: 018 612 1910

Switzerland

Tel: 0848 8028 12
Fax: 0848 8028 13

UK

Tel: 0800 616928
Fax: 0800 616927

USA

Tel: +1 800 526 3593
Fax: +1 877 295 8102

<http://www.gehealthcare.com/lifesciences>

GE Healthcare UK Limited
Amersham Place, Little Chalfont,
Buckinghamshire,
HP7 9NA UK



imagination at work