

U
RANGE



Ultra-Low Freezers

B Medical Systems | U Range

Ultra-low freezers are devices intended for the storage of blood components and blood plasma, human cells, tissues, live virus vaccines, and other laboratory samples at ultra-low temperatures down to -86°C . The devices include an integrated alarm system that warns against unexpected temperature excursions and power failures.

Compliant to DIN 13277 | Medical Device according to MDR (EU) 2017/745, Class IIa or 21CFR Part 864.9700, Class II



medical
systems

SAVING LIVES
THROUGH RELIABLE
AND INNOVATIVE
TECHNOLOGY

Blood Management Solutions



Safety Standards

B Medical Systems

The **Safety Standards developed by B Medical Systems** define certain significant technical features of a product. These ensure the safe storage of the preparations as well as setting the highest standards of safety for the user.

PREMIUM LINE	
●	B Medical Systems Electronics
●	B Medical Systems 7" full touchscreen display
●	Power indicator light and digital temperature indicator (display: 0.1 digits)
●	Self-contained alarm system with integrated battery takes over the alarm function and temperature value measurements in case of system failure for at least 48 hours
●	Acoustic and visual alarm signal in case of temperature alarm and system failure
●	The alarm history function on the electronic stores all the relevant values during a temperature alarm, such as: min., max., average temperature and alarm duration
●	Remote transmission alarm signal (via potential-free contact) in case of temperature alarm (change-over contact)
●	Door opening alarm (visual / acoustic)
●	Designed and tested for climatic class SN / T (ambient temperature range +10°C to +43°C)
●	Interior made from stainless steel
●	Additional remote transmission alarm signal (via potential-free contact) in case of system failure (change-over contact)
●	Smooth castors for optimum flexibility of movement
○	External water cooling
○	Ambient temperature sensor
●	Ethernet interface for the visualization of all operating and control functions (hardware and software settings) via *B Connected monitoring software on a peripheral device (computer)
○	*B Connected - Universal software for the monitoring of refrigeration devices, including the acquisition, recording and visualization of temperature data
●	Integrated remote temperature monitoring device (RTMD), offering real-time worldwide remote monitoring, data access over WEB and GPS position

Ultra-Low Freezers

B Medical Systems | U Range

5 models • Volume 232 > 949 L • Set temperature -82°C • Climate class SN/T • Compliant to DIN 13277 | MDR (EU) 2017/745, Class IIa

In conformity with national and international guidelines, regulations for Medical Devices offering reliability, efficiency and safety at an optimal price.

MODELS U201-401-501-701-901 | PREMIUM LINE



Exclusive integrated electronics handle bar

All functionalities are easily accessible:

- Digital control with full functionalities - B Medical Systems 7" full touchscreen display integrated at optimal level in the door handle with pre-installed connection allowing exclusive °B Connected monitoring functionalities
- Open/close - Heavy duty door lock mechanism



Versatile modularity

Excellent storage capacity and modularity - Large choice of racks offering modularity for every need.



Special features

- Integrated pressure release valve in the door permitting faster opening in one move, and offering easy access for cleaning
- Insulated inner doors for significantly lower loss of cooling when open
- Full compatibility of inner modules and rails, no need to buy new racks, you can keep racks previously used in different models



High quality materials

High quality coating, certified medical devices quality and antibacterial, high quality steel for better longevity and easy hygiene control.



Ultra-low freezers are devices intended for the storage of blood components and blood plasma, human cells, tissues, live virus vaccines, and other laboratory samples at ultra-low temperatures down to -86°C . The devices include an integrated alarm system that warns against unexpected temperature excursions and power failures. U models reflect the highest and most uncompromising requirements with state-of-the-art technology and economy. The cooling system is optimally designed with respect to energy consumption as well as to the development of waste heat and noise. Standard factory setpoint is -82°C , temperature is adjustable from -86°C to -20°C .



Technical Data

General features



U201

U401

U501

U701


U901

<i>Gross / Net volume (l)</i>	232 / 217	478 / 454	634 / 602	791 / 751	949 / 900	
<i>Storage capacity</i>	<i>Vials</i>	15000 (2ml)	34500 (2ml)	46000 (2ml)	57500 (2ml)	69000 (2ml)
	<i>Cryoboxes</i>	150 (model H50)	345 (model H50)	460 (model H50)	575 (model H50)	690 (model H50)
	<i>Plasma bags</i>	150 (350ml)	300 (350ml)	400 (350ml)	500 (350ml)	600 (350ml)
<i>Set temperature (preset)</i>	-82°C					
<i>Set temperature (setting range)</i> <i>can be adjusted in steps of 0.1°C</i>	-86°C to -20°C					
<i>Preset cold / warm alarm limit</i>	-87°C / -77°C					
<i>Climate class (ambient temperature range)</i>	SN / T (+10°C to +43°C)					
<i>Defrosting technique</i>	Manual					
<i>Refrigerant type</i>	R290 / R170					
<i>External dimensions H x W x D (mm)</i>	1293 x 699 x 1039	1988 x 699 x 1039	1988 x 845 x 1039	1988 x 992 x 1039	1988 x 1139 x 1039	
<i>Inner dimensions H x W x D (mm)</i>	657 x 447 x 738	1375 x 447 x 738	1375 x 593 x 738	1375 x 740 x 738	1375 x 887 x 738	
<i>Net weight with standard equipment (kg)</i>	195	250	276	297	320	

Energy Saver	<i>Supply voltage (V)</i>	230	220	230	220	230	220	230	220		
	<i>Frequency (Hz)</i>	50	60	50	60	50	60	50	60		
	<i>Power (W)</i>	2000	2000	2000	2000	2000	2000	2000	2000		
	<i>Energy consumption (kWh/24h)</i>	13.5	13.5	11.0	12.5	11.3	12.5	11.9	13.5	12.4	15.5
	<i>Heat emission (Kcal/h)</i>	209	338	326	405	404	433	426	511	444	567
	<i>Compressor running time (%)</i>	48	48	46	44	51	47	50	55	59	61
	<i>Noise level (dB(A))</i> <i>(at 1m height & 1m distance)</i>	49	50	49	50	52	53	48	49	51	52
<i>Hold over time (-80°C to -60°C)</i>	1 h	1 h	1 h 18 min.	1 h 12 min.	1 h 18 min.	1 h 24 min.	1 h 18 min.	1 h 30 min.	1 h 24 min.	1 h 42 min.	



Technical Data Specifications

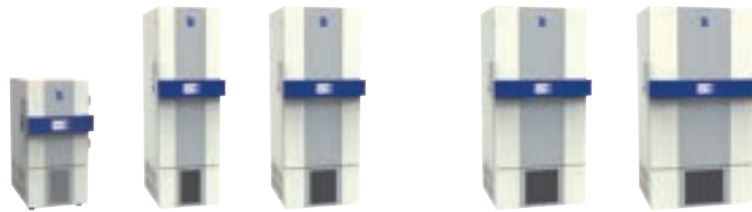

U201-401-501

U701-901

<i>Regulation sensor</i>	PT1000, at 0°C ± 0.1°C
<i>Display sensor in reference body / 100 ml DOW cornring 200-5CST (Silicon Oil)</i>	PT1000, at 0°C ± 0.1°C
<i>Accu data / function time of the control panel when system failure</i>	12 V - 10.0 AH / 48h
<i>Relative humidity at +32°C</i>	≤ 75%
<i>Door insulation (polyurethane)</i>	65-82 mm PU with 20 mm VIP
<i>Cabinet insulation (polyurethane)</i>	77-82 mm PU with 20-30 mm VIP
<i>Material</i>	
<i>Inner body / door</i>	Stainless steel (1.4301)
<i>Outer body / door</i>	DX51D + Z100 coated
<i>Separate interior doors</i>	Anodized aluminium
<i>European Medical Device Regulation</i>	MDR (EU) 2017/745, Class IIa (models: 230 V - 50 Hz)
<i>EMC directive</i>	2014 / 30 / EU (models: 230 V - 50 Hz)
<i>Low voltage directive</i>	2014 / 35 / EU (models: 230 V - 50 Hz)
<i>FDA regulation Medical Device</i>	21CFR Part 864.9700, Class II (models: 220 V - 60 Hz)
<i>cCSAus certificate</i>	Yes (models: 220 V - 60 Hz)
<i>Energy Star certificate</i>	- Yes (models: 220 V - 60 Hz)
<i>Clean room classification ISO 14644-1</i>	Class 5

EXAMPLE OF "STABLE RUNNING (TEMPERATURE STABILITY)"

For all models, temperature distribution graphs for the stages "Cool Down (Cooling)", "Stable Running (temperature stability)" and "Hold Over (temperature graphs, for example in case of system failure)" are available upon inquiry for various voltages. Determined at 5 measuring points per level or compartment (empty).

U401

Voltage: 230V / 50Hz


Set point: -82°C

Ambient temperature: +25°C

Cold alarm: -87°C

Relative humidity: 45%

Warm alarm: -77°C

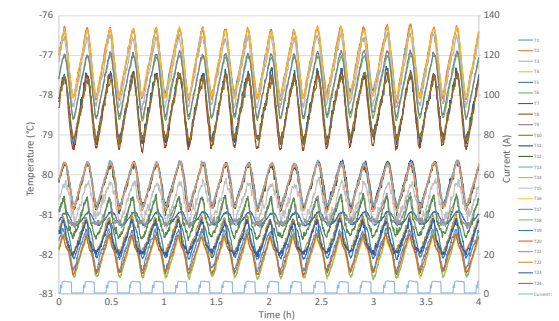



7" FULL TOUCHSCREEN DISPLAY




ERGONOMIC HANDLE



SEPARATE INTERIOR DOORS

Also available with 4 separate interior doors

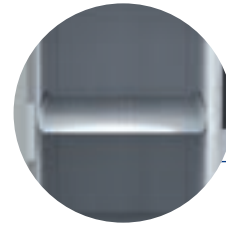

ST-SHELF (STAINLESS STEEL)

Equipment General

	U201	U401	U501	U701	U901
<i>B Medical Systems Electronics</i> 7" full touchscreen display	●	●	●	●	●
<i>Ergonomic handle</i> open/close - heavy duty door lock mechanism	●	●	●	●	●
<i>Integrated pressure release valve</i>	●	●	●	●	●
<i>Separate interior doors</i> in order to minimize the loss of cold air	● 1	● 2 ○ 4	● 2 ○ 4	● 2 ○ 4	● 2 ○ 4
<i>ST-Shelf</i>	● 1	● 3	● 3	● 3	● 3
<i>Ethernet interface</i>	●	●	●	●	●
<i>°B Connected - Monitoring Software</i>	○	○	○	○	○
<i>Integrated remote temperature monitoring device (RTMD)</i>	●	●	●	●	●
<i>Ambient / condenser temperature sensor</i>	○	○	○	○	○
<i>Potential-free alarm contact</i> in case of system failure	●	●	●	●	●
<i>Integrated port for external sensor</i> (installed by customer)	●	●	●	●	●
<i>External water cooling</i>	○	○	○	○	○
<i>CO2 BackUp / Emergency cooling system</i>	○	○	○	○	○
<i>Smooth castors with stabilizers</i>	● 4	● 4	● 4	● 4	● 4
<i>Door hinge Right Left</i>	● ○	● ○	● ○	● ○	● ○



INTEGRATED PRESSURE RELEASE VALVE

This aluminium die cast handle offers a strong and ergonomic feeling that allows the user to open the door with one hand. Inside the freezers, the vacuum will be automatically removed through a release valve, located in the door, in order to allow easy cleaning.

The handle is equipped with an electric lock that can be deactivated through the user interface.

● Standard ○ Optional - Not available


Subject to change without prior notice.




Equipment Accessories




CRYO-RACK
Available For Cryoboxes H50 / H75 / H100

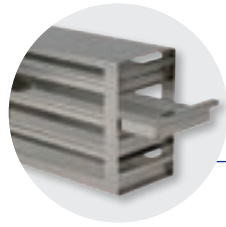

D-RACK
Available For Cryoboxes H50 / H75 / H100

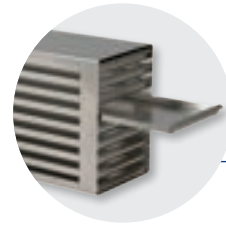

D-RACK MICROTITER

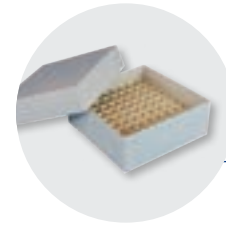


CARTON CRYOBOX
(Shown with carton separator)

	U201	U401	U501	U701	U901	
<i>Large Cryo-Rack H50 for the storage of 30 Cryoboxes H50</i>	-	<input type="radio"/> 9 (max. per unit)	<input type="radio"/> 12 (max. per unit)	<input type="radio"/> 15 (max. per unit)	<input type="radio"/> 18 (max. per unit)	
<i>Cryo-Rack H50 / H75 / H100 for the storage of 25 Cryoboxes H50 / 15x H75 / 10x H100</i>	<input type="radio"/> 6 (max. per unit)	<input type="radio"/> 12 (max. per unit)	<input type="radio"/> 16 (max. per unit)	<input type="radio"/> 20 (max. per unit)	<input type="radio"/> 24 (max. per unit)	
<i>D-Rack H50 / H75 / H100 for the storage of 20 Cryoboxes H50 / 15x H75 / 10x H100</i>	<input type="radio"/> 6 (max. per unit)	<input type="radio"/> 12 (max. per unit)	<input type="radio"/> 16 (max. per unit)	<input type="radio"/> 20 (max. per unit)	<input type="radio"/> 24 (max. per unit)	
<i>D-Rack microtiter H22</i>	<input type="radio"/> 6 (max. per unit)	<input type="radio"/> 12 (max. per unit)	<input type="radio"/> 16 (max. per unit)	<input type="radio"/> 20 (max. per unit)	<input type="radio"/> 24 (max. per unit)	
<i>Carton Cryobox</i>	-	<input type="radio"/> 345 (with 9 Large Cryo-Racks and 3 Cryo-Racks)	<input type="radio"/> 460 (with 12 Large Cryo-Racks and 4 Cryo-Racks)	<input type="radio"/> 575 (with 15 Large Cryo-Racks and 5 Cryo-Racks)	<input type="radio"/> 690 (with 18 Large Cryo-Racks and 6 Cryo-Racks)	
	<i>H50 (H 50 x W 136 x D 136 mm)</i>	<input type="radio"/> 150 (with 6 Cryo-Racks)	<input type="radio"/> 300 (with 12 Cryo-Racks)	<input type="radio"/> 400 (with 16 Cryo-Racks)	<input type="radio"/> 500 (with 20 Cryo-Racks)	<input type="radio"/> 600 (with 24 Cryo-Racks)
	<i>H75 (H 75 x W 136 x D 136 mm)</i>	<input type="radio"/> 120 (with 6 D-Racks)	<input type="radio"/> 240 (with 12 D-Racks)	<input type="radio"/> 320 (with 16 D-Racks)	<input type="radio"/> 400 (with 20 D-Racks)	<input type="radio"/> 480 (with 24 D-Racks)
<i>H100 (H 100 x W 136 x D 136 mm)</i>	<input type="radio"/> 90 (with 6 Cryo-/D-Racks)	<input type="radio"/> 180 (with 12 Cryo-/D-Racks)	<input type="radio"/> 240 (with 16 Cryo-/D-Racks)	<input type="radio"/> 300 (with 20 Cryo-/D-Racks)	<input type="radio"/> 360 (with 24 Cryo-/D-Racks)	
<i>Carton separator for Cryoboxes H50 / H75 / H100</i>	<input type="radio"/> 60 (with 6 Cryo-/D-Racks)	<input type="radio"/> 120 (with 12 Cryo-/D-Racks)	<input type="radio"/> 160 (with 16 Cryo-/D-Racks)	<input type="radio"/> 200 (with 20 Cryo-/D-Racks)	<input type="radio"/> 240 (with 24 Cryo-/D-Racks)	
<i>S-Rack (H 100 x W 113 x D 638 mm) for the storage of bags</i>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
<i>Separator for S-Rack to facilitate the handling and storage management of bags</i>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
<i>Temperature chart recorder pre mounted (in the form of a circular chart recorder)</i>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	




VERSATILE MODULARITY

Excellent storage capacity and modularity - Large choice of racks offering modularity for every need.

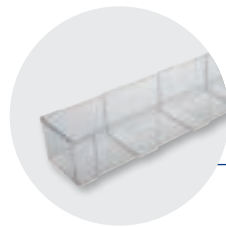

S-RACK

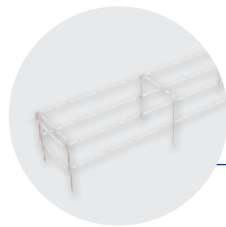

SEPARATOR FOR S-RACK



TEMPERATURE CHART RECORDER
(sample image)

Full functionalities at a glance

B Medical Systems | Electronics

The B Medical Systems Electronics

offers a wide range of adjustment and diagnostic functionalities as well as additional protection / warning operations (via external alarm operations, histories and individual display signals).

PREMIUM LINE | Multifunction electronics with 7" full touchscreen display, integrated at optimal level in the door handle with pre-installed connection allowing exclusive °B Connected monitoring functionalities.



Power switch OFF / ON

- To switch on press the ON/OFF button for 1 second. The green LED will light up
- Switch off protected by password. The door will automatically unlock before powering down



Touchscreen 7" full display

During normal operation the display shows:


- At the top: current voltage, current time, gasket heater and login status (user level)
- Left-hand side: temperature of the regulation sensor (displaying this temperature can be disabled through the settings) and internal temperature displayed in °C (Celsius) or °F (Fahrenheit) in increments of 0.1°
- Right-hand side: touch buttons to access door unlock, clean, menu, alarm mute and message functions



USB port & SD card slot

Acoustic and visual alarm signal

When an alarm occurs, the red alarm LED lights up and an acoustic signal sounds. The corresponding alarm message appears on the display together with the inside measured temperature. If several alarms occur simultaneously, the messages will be shown alternately.

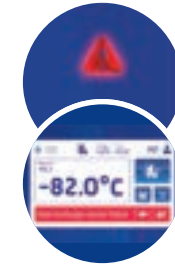

Keycard reader

Login / Unlock using an NFC card

Touch buttons

- DOOR UNLOCK button: unlocks the door of the device
- CLEAN button: the cleaning button allows the user to disable the touch capability of the screen for a short period of time in order to clean the screen. This button is replaced by the logout button when user is logged in
- MENU button: allows access to the menu screen
- ALARM MUTE button: deactivates the acoustic alarm for a predefined period of time
- ACTIVE ALARMS AND WARNINGS button: during a warning or alarm situation, the color of this button and the corresponding level will change. Information on an active alarm and warning is accessible by pressing the message button



Optimum control and protection

B Medical Systems | °B Connected & RTMD



°B CONNECTED | MONITORING SOFTWARE

Universal software for the monitoring of refrigeration devices, including the acquisition, recording and visualization of temperature data.



- Unique monitoring software for the full range of Blood Management Solutions and Medical Refrigeration products
- Web-based interfaces for computers and mobile devices
- Modern design for simple and intuitive use
- Graphic display of temperature curves
- Integrated event and activity history of appliances' components
- Data recording on centralized database for long-term archiving
- Easy setting of specific alarm, via email or SMS alerts
- Generation of reports compiling crucial data and activities
- Temperature and detailed device data export for third-party software
- Important cost advantage compared to a traditional circular chart recorder and its spare parts
- REST API to access raw data directly from the database in read-only mode

KEY BENEFITS:

- Unique interface for the monitoring of the full range of refrigeration devices
- Centralized database providing data access to entire customer network

1 CONFIGURATION




USER MANAGEMENT

- Password login
- Access right



DEVICE SETTINGS

- Device information
- Device & user configuration

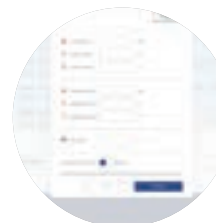

ALARM PROFILE

- Parameters of the alarm profile (day, hour, minute)
- Alarm management

2 DISPLAYING



DETAILED DISPLAY


- Optimized design for mobile and touch compatible devices
- Support of multiple languages
- Real-time data: Current temperatures, set point, upper and lower limit, door state, compressor runtime
- Log data: Events, alarms within date/time, event type, description of the event. Event and alarm confirmation function
 - Parameter settings of each device: Factory settings, real-time clock, language, firmware version, RS485 address, IP address

- °B Connected & RTMD offer a highly flexible system that is adaptable to specific customer requirements
- > Complete & legally safe documentation of temperature data
- > Comprehensive applications and diagnostic functionalities

RTMD | REMOTE TEMPERATURE MONITORING DEVICE (INTEGRATED)

Hardware module monitoring all operating conditions and passing them through our web based Vaclog network ensuring a high performing and reliable cold chain 24/7.

- Real time monitoring of temperature(s), door / lid openings & GPS position
- Includes a SIM chip with prepaid communication for the duration of the warranty period contractually agreed
- Crucial for the quality surveillance of the cold chain and monitoring of vaccines
- Rechargeable battery
- Worldwide remote monitoring & data access over WEB
- Only 2G or 3G coverage is necessary to operate
- Google Maps positioning using integrated GPS module
- Alarms include temperature deviations and door / lid openings
- Alarms are sent by text messages or emails

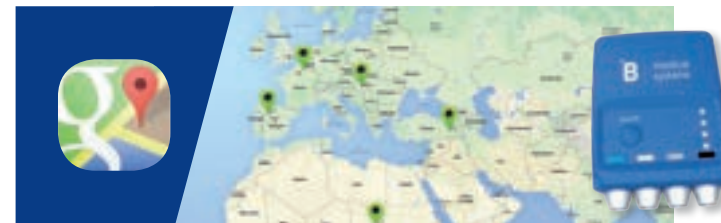


3 ANALYSING




REAL TIME GRAPHICS

- Graphic temperature progress (in real time)
- Adjustable time axis
- Offline navigation

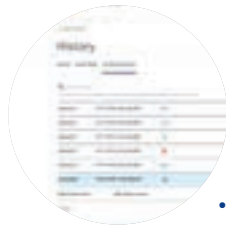


INTERACTIVE ANALYSIS

- No limitation in terms of time range
- Freely scalable
- Detailed information for any point of time
- Temperature data, events and alarms
- Login history and audit trails

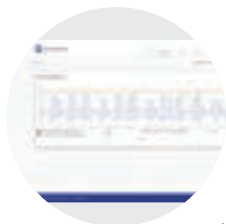

CONTROL CENTER

- Overview (Control Monitor) of all monitored appliances:
- Information regarding status, appliance name, location, connection, upper and lower limits, current temperature, events / alarms
 - Dashboard to combine several temperature measurements in one single graph

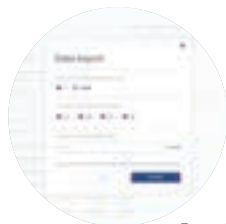

PRINTED ANALYSIS

- Customizable print reports on daily, weekly and monthly basis
- Exports in various file types
- Standardized format



Blood Management Solutions

Safe global blood management: from collection to transfusion, transportation, processing and storage



Vaccine Cold Chain

Reliable solutions for safe vaccination around the world



Medical Refrigeration

State-of-the-art technology for the exacting needs of the medical world



After Sales support and service

We strive to provide you with the highest standards of service; not only through our selected distributors and partners for all your maintenance and service but also our second line trouble shooting and after sales service. This factory-based group of engineers is there to help our partners and yourself to get the best solution for your cold storage needs.

Our Global Expertise



SAVING LIVES THROUGH RELIABLE AND INNOVATIVE TECHNOLOGY

B Medical Systems (formerly Dometic Medical Systems) has more than 40 years' experience in the medical refrigeration sector.

The company, formerly known as Electrolux Medical Systems, was founded in 1979 when the World Health Organization approached Electrolux in Vianden, Luxembourg, to create a solution for the safe storage and transport of vaccines around the world. In 2001, Electrolux Medical Systems became part of the Dometic Group, and was renamed Dometic Medical Systems. Having established a legitimate reputation in the medical equipment industry, the company has also become a global leader in vaccine cold chain.

Luxembourg, in the heart of Europe



B Medical Systems S.à r.l.

17, op der Hei
L - 9809 Hosingen, Luxembourg

Tel.: (+352) 92 07 31-1
Fax: (+352) 92 07 31-300
info@bmedicalsistemas.com



www.bmedicalsistemas.com

WE SUPPORT



Since 2019 B Medical Systems has been committed to the UN Global Compact corporate responsibility initiative and its principles in the areas of human rights, labour, the environment and anti-corruption.