

product codes:

**E70092Y, E70092Z, E70092X**

# Shrimp Alkaline Phosphatase

Removes 5'-phosphates from DNA and RNA


Shrimp Alkaline Phosphatase (SAP) is a Tested User Friendly™ enzyme from USB Corporation. It is a high specific activity, heat-labile alkaline phosphatase that is useful in many applications, including preparing PCR products for labelling, cloning, sequencing (1–3), and SNP analysis. SAP is easily inactivated by heat and offers one-step degradation of dNTPs.

## Dephosphorylation

Alkaline phosphatases can be used for the dephosphorylation of 5'-phosphorylated ends of DNA or RNA for subsequent labelling with <sup>32</sup>P using [ $\gamma$ -<sup>32</sup>P]ATP, T4 polynucleotide kinase, and OptiKinase™. Dephosphorylation also prevents religation of linearized plasmid DNA in cloning experiments. Shrimp Alkaline Phosphatase has approximately the same specific activity as the calf intestine enzyme (800–1000 units/mg at 25 °C, pH 9.6) but, unlike the calf enzyme, SAP can be completely and irreversibly inactivated by heating for 15 min at 65 °C (Fig. 1). No further treatment is necessary.

Shrimp Alkaline Phosphatase is particularly useful in preparing PCR products for downstream applications involving sequencing, SNP analysis, or labelling methods. Typically, the excess dNTPs remaining after PCR interfere with subsequent reactions involving DNA synthesis. Shrimp Alkaline Phosphatase eliminates this problem by dephosphorylating all remaining dNTPs from the PCR mixture in one easy step.

For PCR cleanup, Shrimp Alkaline Phosphatase may be combined with Exonuclease I for removal of residual primers and extraneous single-stranded DNA reaction products. Hence, the use of alternative purification methods, such as columns, gels or magnetic separations, are completely eliminated. For convenience, refer to ExoSAP-IT™, which includes both enzymes in a format that is ready to use.

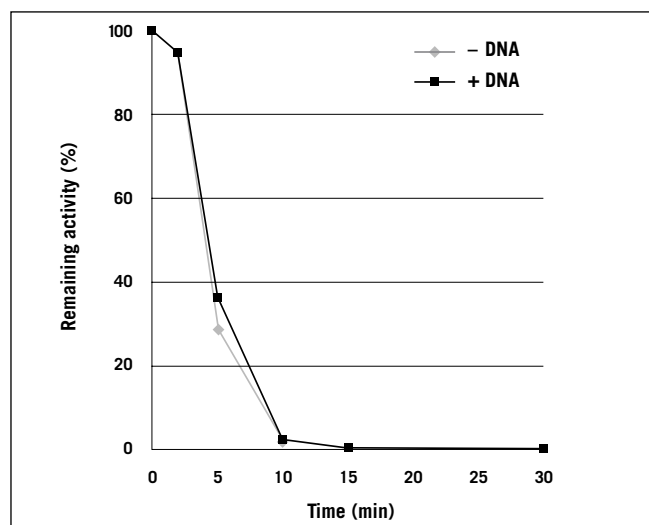

**Fig 1.** SAP inactivation at 65 °C. Reactions were set up using 20 units SAP (excess) with or without 2 µg lambda DNA *EcoR I/Hind III* fragments. The reaction was incubated at 65 °C to inactivate the SAP. Aliquots from the reaction were placed on ice at selected time intervals and assayed for activity in the standard assay.

## Specifications

Shrimp Alkaline Phosphatase is supplied with 10× reaction buffer and SAP dilution buffer; 1 ml each.

**Tested User Friendly Functional Test:** Dephosphorylation of a restriction enzyme digested plasmid (5–20 pmol of 5' ends, 0.1–0.5 units/pmol 5' ends). Reduces re-ligation to < 0.5% compared to the untreated control.

**Purity:** This enzyme is purified to apparent homogeneity and is free of all contaminating endonucleases, exonucleases, and ribonucleases.

**Storage buffer:** 25 mM Tris-HCl (pH 7.5), 1 mM MgCl<sub>2</sub>, 0.1 mM ZnCl<sub>2</sub>, 50% glycerol.

**Concentration:** 1 unit/ul

## Reference

1. Ruan, C. C., *et al. Comments* 17, No. 1, United States Biochemical Corporation (1990).
2. Werle, E. *et al.* Convenient single-step, one tube purification of PCR products for direct sequencing. *Nucleic Acids Res* 22, 4354–4355 (1994).
3. Hanke M. and Wink M. Direct DNA sequencing of PCR-amplified vector inserts following enzymatic degradation of primer and dNTPs. *Biotechniques* 17, 858–860 (1994).

## Ordering information

Shop online at [www.amershambiosciences.com](http://www.amershambiosciences.com)

Shrimp Alkaline Phosphatase	500 units	E70092Y
Shrimp Alkaline Phosphatase	1000 units	E70092Z
Shrimp Alkaline Phosphatase	5000 units	E70092X
Exonuclease I	2500 units	E70073Z
Exonuclease I	5000 units	E70073X
ExoSAP-IT	100 reactions	US78200
ExoSAP-IT	500 reactions	US78201
ExoSAP-IT	2000 reactions	US78202

**Asia Pacific** Tel: +852 2811 8693 Fax: +852 2811 5251 **Australasia** Tel: +61 2 9899 0999 Fax: +61 2 9899 7511 **Austria** Tel: 01 576 0616 19 Fax: 01 576 0616 27  
**Belgium** Tel: 0800 73 888 Fax: 03 272 1637 **Canada** Tel: +1 800 463 5800 Fax: +1 800 567 1008 **Central, East, and Southeast Europe** Tel: +43 1 982 3826  
 Fax: +43 1 985 8327 **Denmark** Tel: 45 16 2400 Fax: 45 16 2424 **Finland & Baltics** Tel: +358 (0)9 512 39 40 Fax: +358 (0)9 512 39 439 **France** Tel: 01 69 35 67 00  
 Fax: 01 69 41 96 77 **Germany** Tel: 0761 4903 490 Fax: 0761 4903 405 **Italy** Tel: 02 27322 1 Fax: 02 27302 212 **Japan** Tel: +81 3 5331 9336 Fax: +81 3 5331 9370  
**Latin America** Tel: +55 11 3933 7300 Fax: +55 11 3933 7304 **Middle East and Africa** Tel: +30 210 96 00 687 Fax: +30 210 96 00 693 **Netherlands** Tel: 0165 580 410  
 Fax: 0165 580 401 **Norway** Tel: 2318 5800 Fax: 2318 6800 **Portugal** Tel: 21 417 70 35 Fax: 21 417 31 84 **Russia & other C.I.S. & N.I.S.** Tel: +7095 232 0250, 956 1137  
 Fax: +7095 230 6377 **Southeast Asia** Tel: +60 3 8024 2080 Fax: +60 3 8024 2090 **Spain** Tel: 93 594 49 50 Fax: 93 594 49 55 **Sweden** Tel: 018 612 1900  
 Fax: 018 612 1910 **Switzerland** Tel: 0848 8028 12 Fax: 0848 8028 13 **UK** Tel: 0800 616928 Fax: 0800 616927 **USA** Tel: +1 800 526 3593 Fax: +1 877 295 8102

### Amersham Biosciences UK Limited

Amersham Place, Little Chalfont, Buckinghamshire, England HP7 9NA

### Amersham Biosciences AB

SE-751 84 Uppsala, Sweden

### Amersham Biosciences Corp.

800 Centennial Avenue, PO Box 1327, Piscataway, NJ 08855 USA

### Amersham Biosciences GmbH

Munzinger Strasse 9, D-79111 Freiburg, Germany

### Amersham Biosciences (SV) Corp.

928 East Arques Avenue, Sunnyvale CA 94085 USA

**Amersham** and **Amersham Biosciences** are trademarks of Amersham plc

**ExoSAP-IT**, **OptiKinase** and **Tested User Friendly** are trademarks of USB Corp.

Exonuclease I/Shrimp Alkaline Phosphatase method of use is covered by one or more of the following patents: 5,756,285 and 5,741,676. ExoSAP-IT is covered by one or more of the following patents: 6,379,940; 6,387,634.

The Polymerase Chain Reaction (PCR) is covered by patents owned by Roche Molecular Systems and F Hoffmann-La Roche Ltd. A license to use the PCR process for certain research and development activities accompanies the purchase of certain reagents from licensed suppliers such as Amersham Biosciences and affiliates when used in conjunction with an authorized thermal cycler.

© Amersham Biosciences Corp. 2003 — All rights reserved

All goods and services are sold subject to the terms and conditions of sale of the company within the Amersham Biosciences group which supplies them. A copy of these terms and conditions is available on request.

[www.amershambiosciences.com](http://www.amershambiosciences.com)

