

Geneseal Instructions for Use

1. Roll Holder Adjustment

(a) Adjust the Geneseal Roll Holder by loosening screws and sliding ends to the width of the roll. Note that film will need to be within the range of the Side Jig Plate (see picture in step 2b).




(b) Insert a roll of Geneseal sealing film. Tighten the screws on the bottom of the roll.




Note: You will probably have to remove the roll to tighten securely.

2. Tray Jig Adjustment

(a) Use these screws to adjust the Genesealer to accommodate your specific tray.

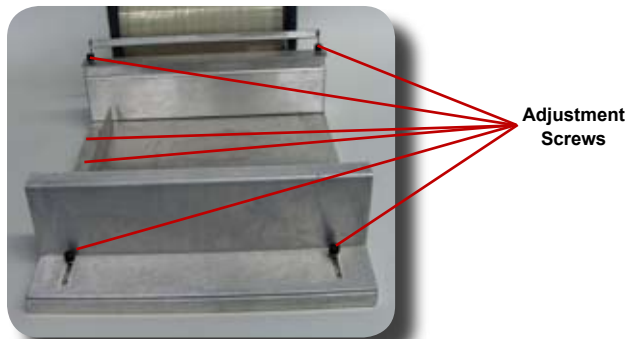

(b) Insert a tray of vials flush with the Front Guide Plate and Side Jig Plate. Note that with assembled corrugated trays, thick end should be oriented forward.

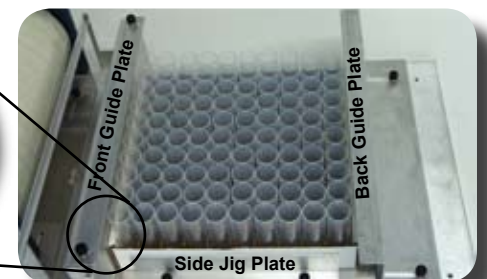


(c) Using the Adjustment Screws, position the adjustable Front and Back Guide Plates such that the edges are nearly touching the vials. This will minimize excess film that results in undesirable overhang on sides of vials.

3. Dispensing Film (after food has cooled ~ 2 hours)

(a) Feed the Geneseal film under the Geneseal Guide and pull to cover the entire tray. Rub film over vials with the flat of the hand to ensure full adhesion. Cut the film with a sharp blade along edges of the Guide Plates. Slide the tray out of the Genesealer. Remove excess film from Back Guide Plate as needed.

Cut along these edges

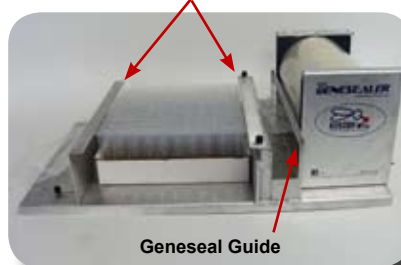

(b) Place another tray (without dividers) or a flat surface (such as a piece of cardboard) onto the sealed vials and flip vials to open-end down and position.

