


Method for Coating Corning® HYPERFlask® Vessels with Collagen or Gelatin

Protocol

CORNING

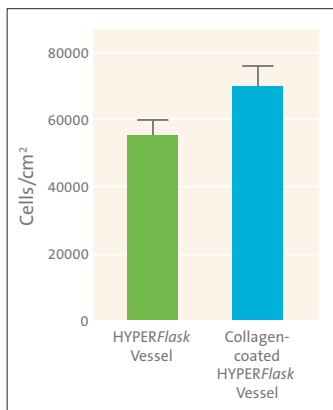

HEK293 cell yields after 4 days of growth on coated and uncoated Corning HYPERFlask vessels. Error bars are the result of 3 separate studies (n = 9).

Introduction

There are many procedures that can be used to coat flasks with collagen, gelatin, or other biological coatings. The following protocols are designed to produce thin coatings to aid in cell attachment and spreading. We recommend using a higher volume of diluted coating solution to make sure the HYPERFlask is completely coated. The final $\mu\text{g}/\text{cm}^2$ will be the same as a standard vessel but a larger volume of more dilute coating solution will be required.

Collagen

Materials

- ▶ Collagen I, rat tail (Corning Cat. No. 354236)
- ▶ Phosphate buffered saline (PBS) (Corning Cat. No. 21-030-CM)
- ▶ Corning HYPERFlask vessel (Corning Cat. No. 10034)

Procedure

NOTE: Due to the viscosity of most collagen, it is recommended to make a diluted stock solution first then use this solution (listed below) to make a working solution at the final desired concentration.

NOTE: Different cell lines will require different collagen coating densities to obtain the desired results. Optimizing the collagen concentration and coating time for your cell line and experimental needs are recommended.

1. Dilute collagen solution to desired concentration with PBS.

Example of calculations to coat one HYPERFlask vessel at $10 \mu\text{g}/\text{cm}^2$ from a 2 mg/mL working solution:

Convert desired coating density (collagen/ cm^2) to collagen concentration (collagen/mL):

$$\frac{1720 \text{ cm}^2 \times 10 \mu\text{g}/\text{cm}^2}{560 \text{ mL}} = 30.71 \mu\text{g}/\text{mL}$$

Calculate total collagen needed to make 560 mL of coating solution:

$$560 \text{ mL} \times 30.71 \mu\text{g}/\text{mL} = 17200 \mu\text{g of collagen}$$

Calculate volume of stock solution needed:

$$\frac{17.2 \text{ mg}}{2 \text{ mg}/\text{mL}} = 8.6 \text{ mL of 2 mg}/\text{mL stock solution}$$

Combine: 8.6 mL stock solution with 551.4 mL PBS = 560 mL.

2. Add 560 mL of diluted collagen solution to the HYPERFlask vessel to fill it completely.

NOTE: To achieve a uniform coating, the HYPERFlask vessel must be filled completely.

3. Incubate the HYPERFlask vessel at room temperature for the desired length of time (typically 2 hours).

NOTE: If a 37°C incubation is preferred, the HYPERFlask vessel must be filled completely to prevent pressure buildup within the vessel from the air heating and expanding.

4. After incubation, aspirate or pour collagen solution from the HYPERFlask vessel and wash once with 50 mL of PBS, making sure that all layers of the HYPERFlask are washed appropriately. The vessel is now ready for use or storage at 4°C until needed.

Gelatin

Materials

- ▶ Porcine gelatin (MilliporeSigma Cat. No. G1890-100G)
- ▶ Sterile distilled deionized water (Corning Cat. No. 25-055-CV)
- ▶ Corning® HYPERFlask® vessel (Corning Cat. No. 10034)

Procedure

NOTE: Different cell lines may require different gelatin coating concentrations to obtain the desired results. Optimizing the concentration and coating time for your cell line and experimental needs is recommended.

1. Make gelatin solution to desired concentration with water.

Example of calculations to coat one HYPERFlask vessel with a 0.1% gelatin solution (1 g/1L):

Calculate total gelatin needed to make 560 mL of coating solution:

$$560 \text{ mL} \times 1 \text{ g}/1000 \text{ mL} = 0.56 \text{ g of gelatin}$$

2. Add gelatin to the water. Gently heating the solution will help the gelatin to dissolve.
3. Autoclave the gelatin solution to sterilize.
4. Cool the gelatin to room temperature before use or store at 4°C until ready to use.
5. Add 560 mL of room temperature gelatin solution to the HYPERFlask vessel to fill it completely.
6. Incubate the HYPERFlask vessel at 37°C for at least 30 minutes.

NOTE: To achieve a uniform coating and to prevent pressure buildup within the vessel, the HYPERFlask vessel must be filled completely during the 30-minute or longer incubation.

7. After incubation, aspirate or pour out the gelatin solution from the HYPERFlask vessel and use.

For more specific information on claims, visit the Certificates page at www.corning.com/lifesciences.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

For additional product or technical information, email ScientificSupport@corning.com, visit www.corning.com/lifesciences, or call 800.492.1110. Outside the United States, call +1.978.442.2200 or contact your local Corning sales office.

Corning Incorporated
Life Sciences
836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476
www.corning.com/lifesciences

**Worldwide
Support Offices**

ASIA/PACIFIC
Australia/New Zealand
t 61 427286832

China
t 86 21 3338 4338
f 86 21 3338 4300

India
t 91 124 4604000
f 91 124 4604099

Japan
t 81 3-3586 1996
f 81 3-3586 1291

Korea
t 82 2-796-9500
f 82 2-796-9300

Singapore
t 65 6572-9740
f 65 6861-2913

Taiwan
t 886 2-2716-0338
f 886 2-2516-7500

EUROPE

France
t 0800 916 882
f 0800 918 636

Germany
t 0800 101 1153
f 0800 101 2427

The Netherlands
t 020 655 79 28
f 020 659 76 73

United Kingdom
t 0800 376 8660
f 0800 279 1117

**All Other European
Countries**
t 31 (0) 20 659 60 51
f 31 (0) 20 659 76 73

LATIN AMERICA
grupoLA@corning.com

Brasil
t (55-11) 3089-7400

Mexico
t (52-81) 8158-8400

CORNING | **FALCON** | **AXYGEN** | **GOSELIN** | **PYREX**