

# MyPipetman®

SELECT

## IMPROVED PIPETTING AND PERSONALIZATION OPTION

Air-Displacement Pipettes



**Ergonomic, safe, and reliable, MyPIPETMAN® Select offers volume setting choice and optional handle personalization for optimal use in various protocols.**

### Select, Set, and Secure Volumes with the Trilock™ System

The unique and patented Trilock™ volume-locking system enables you to choose different positions when setting the volume:

- Free volume setting position for quick and effortless volume adjustment
- Light ratchet volume setting for precise volume adjustment
- Volume locked position for secure pipetting

With this new Trilock system, the volume is adjusted simply by turning the push button of the pipette.



### Ergonomic Design

Lightweight and well-balanced, MyPIPETMAN Select has an elegant design with a textured effect for a better grip.

With no compromise on quality and performance, the built-in ejection system provides reduced tip ejection forces. This patented ejection system reinforces the pipette's ergonomics for convenient use throughout the day.

### Handle Printing with Your Text on Demand

MyPIPETMAN Select can have their handles printed on demand. Identify your pipette to avoid any risk of contamination by dedicating it to a specific liquid or workflow step and improve your lab Lean Management process.





## Extra Personalization and Accessories

Other accessories like colored windows, tip holders, and tip ejectors are available in different colors and material for additional customization and ease of use. Please visit [www.gilson.com](http://www.gilson.com) for ordering information.



Colored Windows	
Model	Part Number
Window 3digit Dark Blue MyPIPETMAN	FP070101
Window 3digit White MyPIPETMAN	FP070102
Window 3digit Blue MyPIPETMAN	FP070103
Window 3digit Red MyPIPETMAN	FP070104
Window 3digit Yellow MyPIPETMAN	FP070105
Window 3digit Green MyPIPETMAN	FP070106

## Pipette Stands

Different pipette stands and holders allow you to store pipettes vertically, protecting them from contamination.


CARROUSEL™



Universal Multichannel Stand



TRIO™

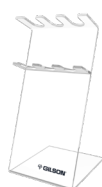


SINGLE™



Accessories	
Description	Part Number
CARROUSEL™ pipette Stand (7 pipettes)	F161401
TRIO™ Pipette Stand (3 pipettes)	F161405
Universal Multichannel Stand	F161417
SINGLE™ Pipette Holder	F161406

## Pipette Kit

MyPIPETMAN Select single channel models are also available in a **4-pipette kit**, including tips and accessories to start an experiment right away.

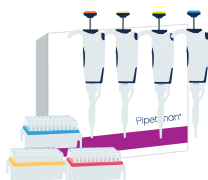


MyPIPETMAN® Select Kit				
Model	Part Number	Models	PIPETMAN® DIAMOND Tips	Accessories
MyPIPETMAN® Select 4-Pipette Kit	F167390	P2 P20 P200 P1000	DL10 D200 D1000	4 SINGLE pipette holders 4 Plastic ejectors 6 colored windows

MyPIPETMAN® Select				
Model	PIPETMAN® DIAMOND Tips		Part Number	Volume (µL)
<b>Single Channel Models</b>				
P2	D10 DL10	DF10ST DFL10ST	FP10001S	0.2 - 2
P10	D10 DL10	DF10ST DFL10ST	FP10002S	1 - 10
P20	D200	DF30ST	FP10003S	2 - 20
P100	D200	DF100ST	FP10004S	10 - 100
P200	D200 D300	DF200ST DF300ST	FP10005S	20 - 200
P1000	D1000 D1200	DF1000ST DF1200ST	FP10006S	100 - 1000
P5000	D5000		FP10007S	500 - 5000
P10mL	D10mL		FP10008S	1000 - 10000
<b>Multichannel Models</b>				
P8x10	D10 DL10	DF10ST DFL10ST	FP10009S	0.5 - 10
P12x10			FP10010S	
P8x20	DL10 D200	DFL10ST DF30ST	FP10011S	2 - 20
P12x20			FP10012S	
P8x200	D200 D300	DF200ST DF300ST	FP10013S	20 - 200
P12x200			FP10014S	
P8x300	D200 D300	DF200ST DF300ST	FP10015S	20 - 300
P12x300			FP10016S	



MyPIPETMAN is also available for Enterprises offering both functional and product customization options. Benefit from our Gilson brand and quality standards to build your own specific line of pipettes.

Please visit [www.gilson.com](http://www.gilson.com) for more information on MyPIPETMAN Enterprise

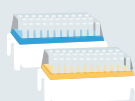

PIPETMAN DIAMOND® Tips are designed to fit your PIPETMAN perfectly and guarantee maximum accuracy and precision.



Benefit from a special care with our flexible, tiered service level offerings, tailored to fit your lab's needs, while increasing the lifespan of your pipettes.

Please visit [www.gilson.com](http://www.gilson.com) for more information or ordering

Please contact your Gilson Account Manager for special text printing. Personalized pipettes will have the same part number with SP at the end.