




Waterbath

WPE 45

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

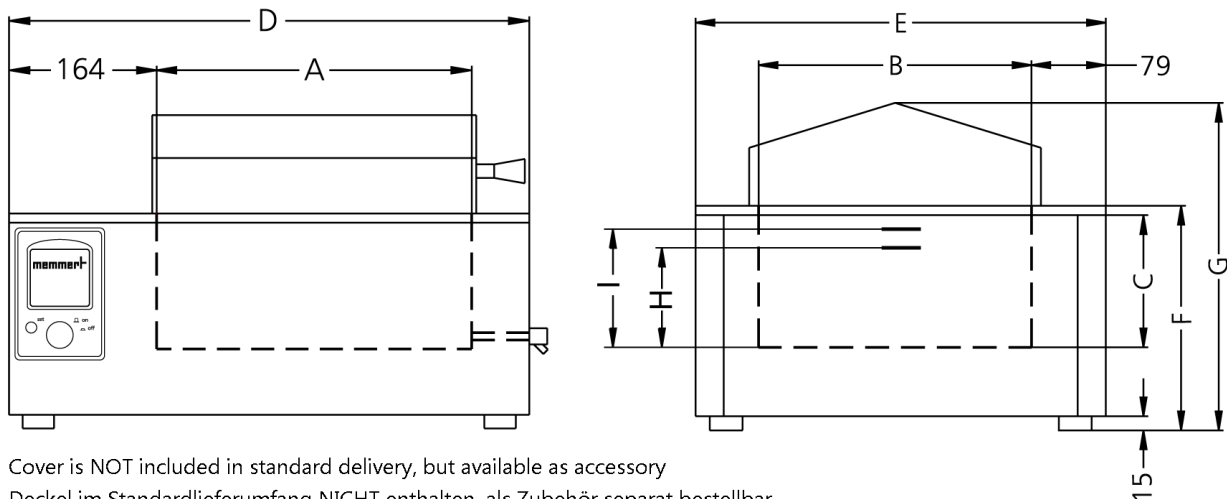

Temperature

Temperature range in °C min. 15°C above ambient up to +95°C

resolution of display and setting accuracy 0,1°C up to +99,9°C / 1°C from 100°C

Control of standard components

Circulation pump pump can be switched on and off

Circulation pump pump is switched off at programme end

Controller digital display (LED) of all set parameters, such as temperature and alarm values (0,1°C resolution) as well as time values

Timer digital timer from 1 min. up to 999 hours for: ON (continuous operation), DELAYED ON, HOLD or HOLD set-temperature dependent with guaranteed holding time

Calibration on controller

Control technology

Controller Electronic microprocessor temperature controller with auto-diagnostic system

Safety

Temperature control 2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value

Temperature control mechanical temperature limiter TB protection class 1 switching the heating off at approx. 30°C above max. temperature of the bath

Temperature control independently working, electronic overtemperature controller TWW protection class 3.1 or overtemperature limiter TWB protection class 2, adjustable in set-up menu by customer

Temperature control in case of overtemperature due to failure, the heating is switched off at approx. 10°C above the set temperature (fixed value)

Autodiagnostic system fuzzy-supported PID microprocessor controller with integrated autodiagnostic system with fault indication

Alarm visual and acoustic alarm at programme end and as input acknowledgement as well as in case of low liquid level, heating is switched off automatically

Heating concept

Heating Baths corrosion-proof large-area heating on three sides

Stainless steel interior

Dimensions W x H x D in mm $w_{(A)} \times h_{(C)} \times d_{(B)}$: 590 x 220 x 350 mm

Interior easy-to-clean interior, made of stainless steel, reinforced by deep drawn ribbing, material no. 1.4301 (ASTM 304), laser-welded

Volume 45 l

Liquid level min. (H) 177 mm

Liquid level max. (l) 200 mm

Textured stainless steel casing

Dimensions $w_{(D)} \times h_{(G)} \times d_{(E)}$: 818 x 401 x 516 mm

Electrical data

Voltage 230 V, 50/60 Hz

Electrical load approx. 2800 W (during heating)

Ambient conditions

Installation The vent openings in the left and back side must remain unobstructed. Minimum wall spacing on all sides is 80 mm. The minimum spacing from the top of the bath to the next ceiling is 750 mm.

Ambient temperature +5 °C to +40 °C

Humidity rh max. 80 %, non-condensing

Overvoltage category II

Pollution degree 2

Packing/shipping data

Transport information The appliances must be transported upright

Customs tariff number 8419 8998

Country of origin Federal Republic of Germany

WEEE-Reg.-No. DE 66812464

Dimensions approx incl. carton w x h x d: 910 x 610 x 460 cm

Net weight approx. 26 kg

Gross weight carton approx. 33 kg

Standard units are safety-approved and bear the test marks

