


**ID-Nr. 30803-050**  
**Model CS.195.105.MM.4WD**

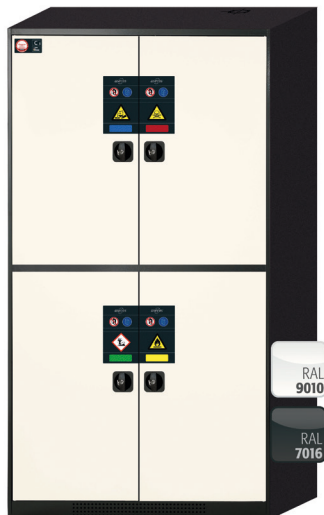

**Cabinet for chemicals CS-CLASSIC-MultiRisk**  
**Model CS.195.105.MM.4WD, body colour anthracite grey**  
**(RAL 7016), with wing doors in pure white (RAL 9010)**

**Safe, approved and separated indoor storage of different non-flammable hazardous substances in one cabinet**

**Function / construction:**

- Robust construction and longevity — sturdy all welded construction suitable for highest industrial demands, high quality, scratch-resistant structured paint finish, double-walled door construction
- **All in one:** different compartments (multiple separation) allow the storage of different hazardous materials in one cabinet
- **No unauthorised use:** doors lockable with cylinder locking (integrated in the turning handle), each compartment can be opened and locked separately
- **Easy alignment:** adjusting aids to compensate for uneven floor
- **Ventilation:** natural ventilation due to ventilation openings at the cabinet bottom, ready for connection to a technical exhaust system

Dimensions	
Width (external)	1055.00 mm
Depth (external)	520.00 mm
Height (external)	1950.00 mm
Width (internal)	475.00 mm
Depth (internal)	495.00 mm
Height (internal)	905.00 mm

The available internal dimensions can vary depending on the interior equipment package.


Weight / Load	
Weight without interior equipment	80 kg
Distributed load	421.00 kg/m <sup>2</sup>
Maximum load	600 kg

Ventilation	
Extraction air	75 DN
Minimum exhaust air quantity 10 times	9 m <sup>3</sup> /h
Differential pressure (10 times air change)	2 Pa


Colouring / Surface finish	
Body colour code	RAL 7016
Body colour	anthracite grey
Body surface finish	textured
Door / Drawer colour code	RAL 9010
Door / Drawer colour	Pure white
Door / Drawer surface finish	textured

General Information	
Number of doors	4
Depth with open doors	1026 mm
Customs tariff number	94032080

Front view

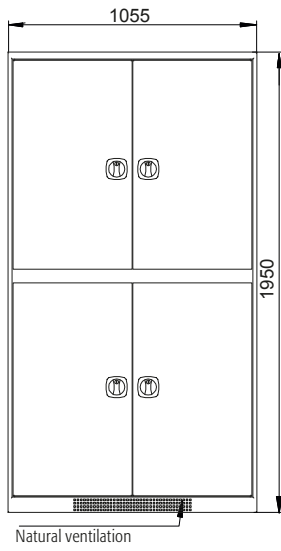

Side view

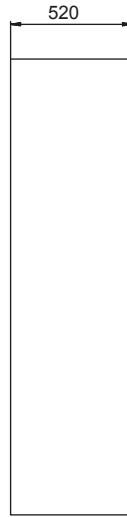

Top view

