

## T309

### CRYOVIAL®

#### External Thread Design with Lip Seal





Made of specially formulated polypropylene

Designed for storing biological material, human or animal cells, at temperatures as low as -196 °C (but should be used only in the gas phase of liquid nitrogen). The cap features a long skirt for one hand aseptic methods, and a super fast thread design that allows tightening or removal with a mere 1¼ turn, and an inside thread design that will not contribute to possible contamination. A specially designed lip inside the cap ensures a leakproof seal even at very low temperatures. Closures and tubes are both made of polypropylene having the same coefficient of expansion, which further enhances the leakproof qualities of these vials at various temperatures. Tubes are provided with a white marking area for sample identification and can be color coded by the use of a CAPINSERT™ (Series T312). T309-2 can be centrifuged up to 17,000g. Vials are sterilized by gamma radiation and are packaged in unique tamperproof, resealable, safety-lock bags of 100. Autoclavable

Tubes on this page are certified RNase, DNase, Pyrogen and DNA-free



For Capinsert™, please refer to T312 on page 98.

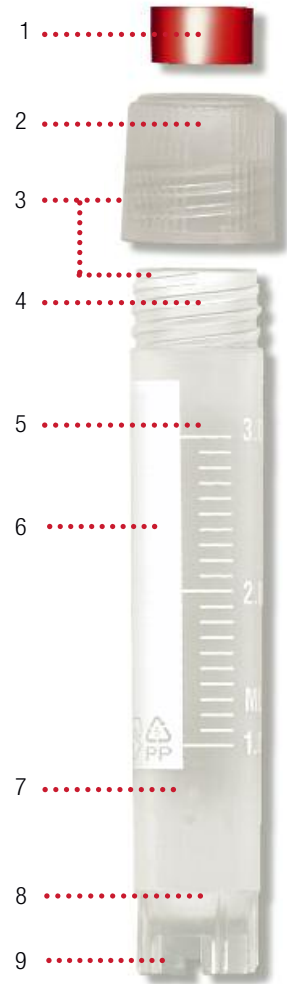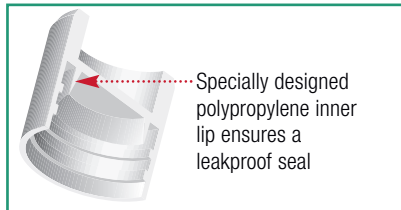

- 1- A CAPINSERT™ is available in 11 different colors / Perfect for color coding (See T312 Series on page 98)
- 2- Vertical ribs facilitate cap removal
- 3- Both cap and tube are made of same polypropylene material, therefore same coefficient of expansion ensures secure seal at all temperatures
- 4- Super fast 1¼ turn thread design
- 5- Thick wall makes vial almost unbreakable
- 6- Large white marking area
- 7- Excellent clarity makes sample easy to see
- 8- Round bottom / Very easy to empty contents completely
- 9- Many sizes available as self-standing with universal locking base

For IVD use

Cat. #	T309-1A	T309-2	T309-2A	T309-3A	T309-4A	T309-5A
Volume (ml)	1.2	2	2	3	4	5
Size (mm)	12.5 x 42	12.5 x 47	12.5 x 49	12.5 x 71	12.5 x 77	12.5 x 91
Self-Standing	•		•	•	•	•
Round Bottom		•				
Qty/Bag	100	100	100	100	100	100
Qty/Cs	1000	1000	1000	1000	1000	1000