

DOSING PUMP FOR LIQUIDS

SIMDOS 02



- ▶ Simple and intuitive operation
- ▶ Easy calibration
- ▶ Cycle mode
- ▶ Adaptable to fluid characteristics
- ▶ Self priming, can run dry
- ▶ Analog and impulse control (RC)
- ▶ Splashproof IP65
- ▶ Small footprint

Fields of application

- Chemistry
- Synthetic organic chemistry
- Pharmacology
- Adhesives and glue
- Polymers & polymer coating
- Food technology research
- Detergents

Hydraulic data

Flow rate	Dosing volume	Suction height	Max. pressure
30 µl/min–20 ml/min	30 µl–999 ml	2 mWg	6 bar

Chemical resistance

Code	Pump head	Diaphragm	Valves	Filter ¹⁾
KT	PP	PTFE covered	FFKM Kalrez [®]	PEEK
TT	PVDF	PTFE covered	FFKM Kalrez [®]	PVDF
FT	PTFE	PTFE covered	FFKM Kalrez [®] *	PEEK
ST	STAINLESS STEEL	PTFE covered	FFKM Kalrez [®] *	PEEK

* On request FFKM Chemraz[®]
 1) material for housing and filter mesh

Functions

S Version	RC Version	RC Plus Version
<ul style="list-style-type: none"> • Manual operation 	<ul style="list-style-type: none"> • Manual operation • Analog control: 0–5 V, 0–10 V, 0–20 mA, 4–20 mA, from 1 to 100% • Start/Stop through logic control (TTL) 	<ul style="list-style-type: none"> • Reset/Prime through logic control (TTL) • Output signal fault alarm • Cable for external control included • Same as RC • Plus RS 232


Technical data

- Main supply 100–240 V / 50–60 Hz
- Ambient temperature allowed +5 to +40 °C
- Liquid temperature allowed +5 to +80 °C
- Maximum viscosity allowed 150 cSt
- Accuracy +/-2% (nominal value)
- Repeatability +/-1%
- Power consumption 10 W
- Protection class IP65
- Dimensions: 150 x 93 x 144
- Weight: 0.9 kg
- Hydraulic connection UNF 1/4"-28

ACCESSORIES

Valves kit

	ID-Nr.
Chemraz® special valves for certain chemicals	164984

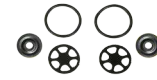

Special Valves CHEMRAZ®

Diaphragm pressure control valve with support when the feeding tank is above the pump

	ID-Nr.
Kit JFDV 30 KTZ	166286
Kit JFDV 30 TTZ	166287
Kit JFDV 30 FTZ	166288



Filter

Typ	Material ¹⁾	Connector	ID-Nr.
FS 60 T	PVDF	UNF 1/4"-28	165210
FS 60 X	PEEK	UNF 1/4"-28	165212

¹⁾ material for housing and filter mesh



Connection tube

	ID-Nr.
Tube 1/16" with connector UNF 1/4"-28, 1 m for fine dosing	166335



Mounting plate

	ID-Nr.
for wall mounting	160473



Standfitting

	ID-Nr.
for retort stand	160474



Footswitch

	ID-Nr.
for impulse start and stop	155872



SIMDOS® is a registered trademark of KNF Flodos AG. All rights reserved.
 DuPont™ Kalrez® is a registered trademark of E. I. du Pont de Nemours and company or one of its subsidiaries. All rights reserved.
 Chemraz® is a registered trademark of Greene, Tweed & Co. All rights reserved.

SALES CENTERS

Germany

KNF Neuberger GmbH
D-79112 Freiburg
Tel. +49 7664 5909 0
info@knf.de
www.knf.de

Benelux Netherlands

KNF Verder B.V.
NL-3451 GG Vleuten
Tel. +31 30 677 92 40
info@knf-verder.nl
www.knf-verder.nl

Benelux Belgium and Luxemburg

KNF Verder N.V.
B-2630 Aartselaar
Tel. +32 3 871 96 24
info@knf.be
www.knf.be

China

KNF Technology (Shanghai) Co., Ltd.
Shanghai 201203
Tel. +86 21 5109 9695
info@knf.com.cn
www.knf.com.cn

France, Morocco, Tunisia, Algeria

KNF Neuberger SAS
F-68128 Village-Neuf
Tel. +33 389 70 35 00
info@knf.fr
www.knf.fr

UK

KNF Neuberger U.K., Ltd.
Witney, Oxfordshire OX28 4FA
Tel. +44 1993 77 83 73
info@knf.co.uk
www.knf.co.uk

India

KNF Pumps + Systems (India) Pvt. Ltd.
Hinjewadi, Pune 411 057
Tel. +91 20 640 13 923
+91 20 640 08 923
info@knfpumps.in
www.knfpumps.in

Italy

KNF ITALIA S.r.l.
I-20063 Cernusco s. Naviglio MI
Tel. +39 02 272 03 860
info@knf.it
www.knf.it

Japan

KNF Japan Co.Ltd.
Tokyo 104-0033
Tel. +81 3 3551 7931
info@knf.co.jp
www.knf.co.jp

Korea

KNF Neuberger Ltd.
135-502, Seoul
Tel. +82 2 959 0255
knf@knfkorea.com
www.knfkorea.com

Sweden, Finland, Denmark, Norway

KNF Neuberger AB
SE-11743 Stockholm
Tel. +46 8 744 51 13
info@knf.se
www.knf.se

Switzerland

KNF Neuberger AG
CH-8362 Balterswil
Tel. +41 71 973 99 30
knf@knf.ch
www.knf.ch

Taiwan

KNF Neuberger Ltd.
Taipei City, 11490
Tel. +886 2 2794 1011
knftwn@knftwn.com.tw
www.knftwn.com.tw

USA, Canada

KNF Neuberger, Inc.
Trenton, New Jersey 08691-1810
Tel. +1 609 890 8600
knfusa@knf.com
www.knfusa.com

Latin America

Tel. +1 609 649 1010
gb@knf.com
www.knf.com/es