



DATA SHEET


## Product Code SLS1244 – SLS Select 2.0ml 96 round well U bottom plate

### About SLS1244

The 2mL round well 'U' bottom plate (SLS1244) has been researched in depth to enhance its fit into the heater blocks featured on automated liquid handlers and a wide range of additional workstations.

The plate is manufactured in SLS Select's class 8 cleanrooms using a medical-grade polypropylene. This polypropylene helps achieve very low concentrations of extractable elements and has a temperature capability range of -80°C to +121°C.

This plate has been developed precisely to conform to the footprint dimensions defined by SLAS/ANSI standards. This guarantees its compatibility with both manual and automated sample handling systems, readers, and microplate washers. The storage plate has stacking features for ease of use in plate and automation hotels.

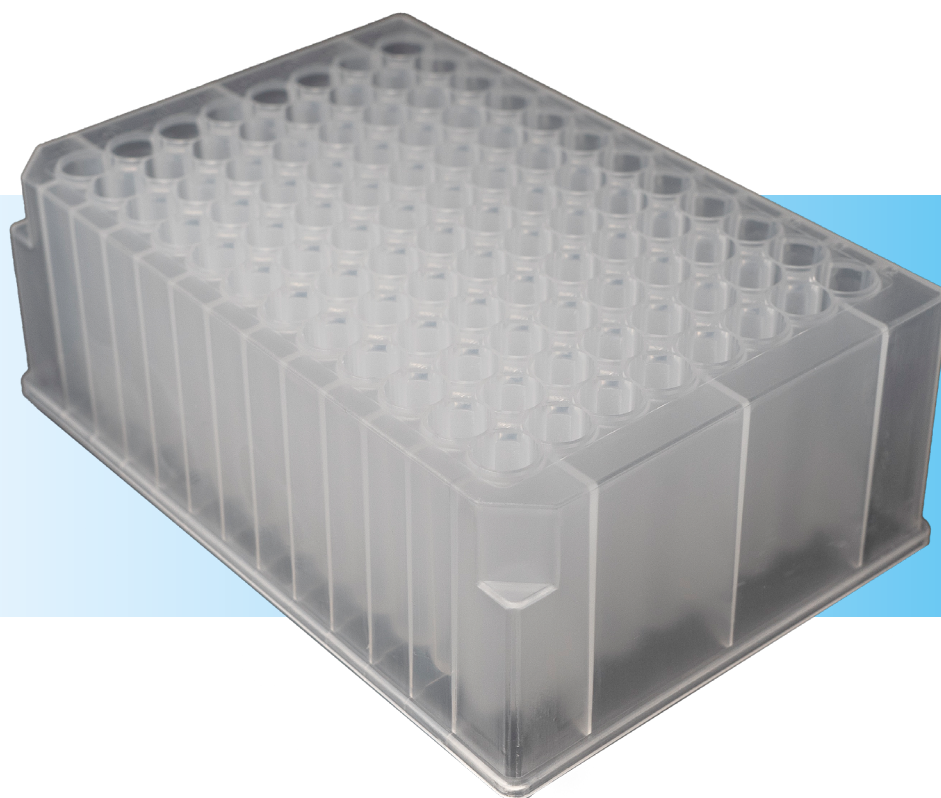

**United Kingdom:**  
slsinfo@scientific-labs.com

**International:**  
slsexport@scientific-labs.com

**Northern Ireland:**  
NI.sales@scientific-labs.com

**Republic of Ireland:**  
ROI.sales@scientific-labs.com

 **SCIENTIFIC  
LABORATORY  
SUPPLIES**



DATA SHEET

## Product Code SLS1244 – SLS Select 2.0ml 96 round well U bottom plate

### KEY FACTS:

- Manufactured in ISO Class 8 Cleanroom
- Alphanumeric well coding
- High temperature and chemical resistance (-80°C to +121°C)
- Autoclavable
- RNase, DNase Free and Non-Pyrogenic
- Stackable for use in plate & automation hotels
- Stack Height for 2 plates 87.90mm
- Available sterile or non-sterile
- Flat surface for sealing with adhesive and heat seals




### ORDERING INFORMATION:

- Case Dimensions: 47 x 47 x 18cm
- Case Weight: 6kg
- Case qty: 50 Plates
- Packsize: 10 Bags x 5 Plates
- Standard pallet quantity: 34 cases

SLS code	Qty per pack	Wells	Well bottom	Max volume mL	L	W	H	Well depth	Well diameter	Stack height of 2 plates	Silicone cap matt
SLS1244	50	96	U	2.0	127.75	85.48	44.5	43	8.5	87.9	SLS1290, SLS1296

SLS1312 available as a sterile version SLS code: SLS1306

### PRODUCT DRAWING: SLS1244

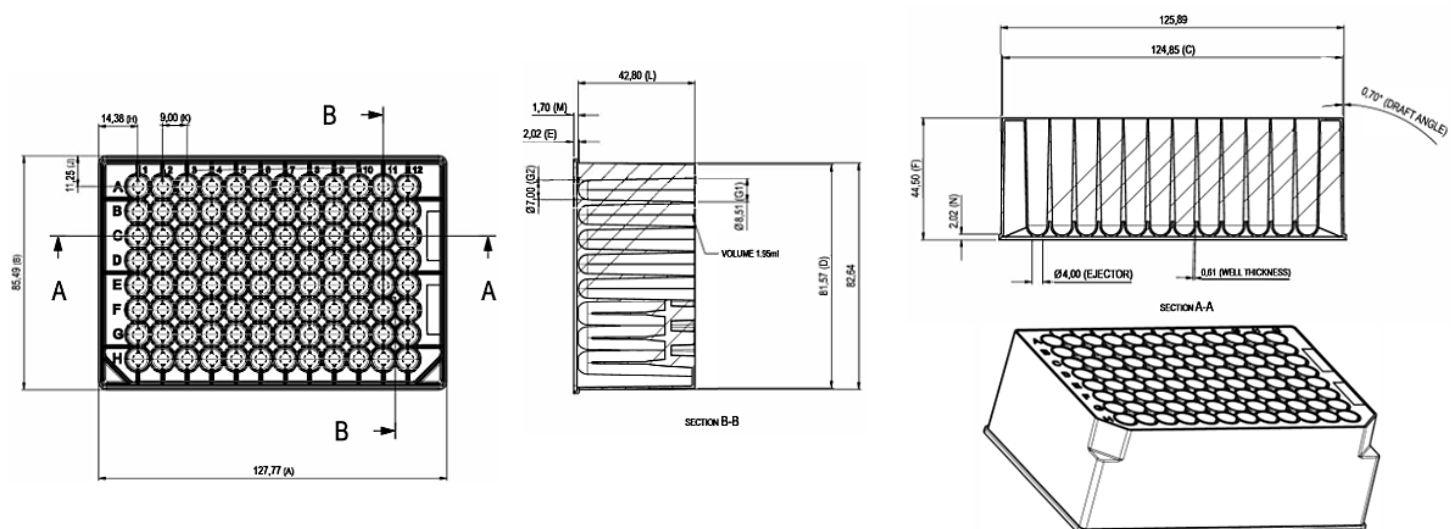

United Kingdom:  
slsinfo@scientific-labs.com

International:  
slsexport@scientific-labs.com

Northern Ireland:  
NI.sales@scientific-labs.com

Republic of Ireland:  
ROI.sales@scientific-labs.com

