



DATA SHEET

Product Code SLS1316 – SLS Select 2.2mL square well ‘U’ bottom plate, with notches

About SLS1316

The 2.2-mL square well ‘U’ bottom plate (SLS1316) has been developed to enable the base of the well to be in contact with heater-shaker blocks which improves the performance of the process.

The plate has the added benefit of notches on the sidewall for compatibility with KingFisher™ platforms.

The plate is manufactured in SLS Select’s class 8 cleanrooms using a medical-grade polypropylene. This polypropylene helps achieve very low concentrations of extractable elements and has a temperature capability range of -80°C to +121°C.

This plate has been developed precisely to conform to the footprint dimensions defined by SLAS/ANSI standards. This guarantees its compatibility with both manual and automated sample handling systems, readers, and microplate washers.

SLS Select also offer a version of this plate without notches. Product code SLS1312.



United Kingdom:
slsinfo@scientific-labs.com

International:
slsexport@scientific-labs.com

Northern Ireland:
NI.sales@scientific-labs.com

Republic of Ireland:
ROI.sales@scientific-labs.com

 **SCIENTIFIC
LABORATORY
SUPPLIES**



DATA SHEET

Product Code SLS1316 – SLS Select 2.2mL square well ‘U’ bottom plate, with notches

KEY FACTS:

- Manufactured in ISO Class 8 Cleanroom
- Alphanumeric well coding
- High temperature and chemical resistance (-80°C to +121°C)
- Autoclavable
- RNase, DNase Free and Non-Pyrogenic
- Stackable for use in plate & automation hotels
- Stack Height for 2 plates 85.81mm
- Available sterile or non-sterile
- Flat surface for sealing with adhesive and heat seals




ORDERING INFORMATION:

- Case Dimensions: 47 x 47 x 18cm
- Case Weight: 6kg
- Case qty: 50 Plates
- Packsize: 10 Bags x 5 Plates
- Standard pallet quantity: 34 cases

SLS code	Qty per pack	Wells	Well bottom	Max volume mL	L	W	H	Well depth	Well diameter	Stack height of 2 plates	Silicone cap matt
SLS1312	50	96	U	2.2	127.75	85.48	44.5	40	8.24	85.81	SLS1284, SLS1286, SLS1294
SLS1316	50	96	U	2.2	127.75	85.48	44.5	40	8.24	85.81	SLS1284, SLS1286, SLS1294

SLS1312 available as a sterile version SLS code: SLS1314
 SLS1316 available as a sterile version SLS code: SLS1318

PRODUCT DRAWING: SLS1316

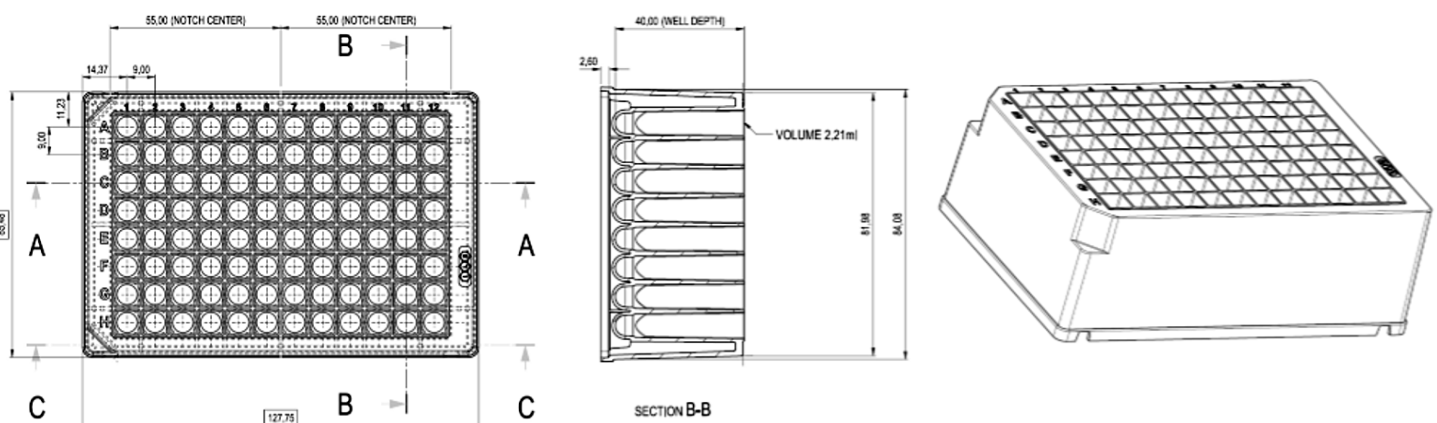

United Kingdom:
 slsinfo@scientific-labs.com

International:
 slsexport@scientific-labs.com

Northern Ireland:
 NI.sales@scientific-labs.com

Republic of Ireland:
 ROI.sales@scientific-labs.com

