

## SLS2046 25Litre Ultrasonic Bath

The SLS Lab Pro range of digital baths offer outstanding performance for use in scientific and laboratory applications, chiropody, veterinary and many other professional environments.

Incorporating Frequency Leap Technology, homogenous cavitation activity is generated throughout the ultrasonic bath reducing dead spots and standing waves. This provides a more consistent scrubbing action for excellent results.

The SLS2046 is a 25Litre heated, digital ultrasonic bath.

The menu driven control panel allows simple adjustment of cycle time (0 to 99 minutes), temperature (ambient to 70°C), ultrasonic power and frequency mode for precise and accurate processing.

Combining features such as degas mode, fluid level sensor and a memory card port, the SLS2046 offers repeatable and traceable cleaning performance with consistent cleaning results every time.




Specifications are subject to change for continual improvements



Made in Britain



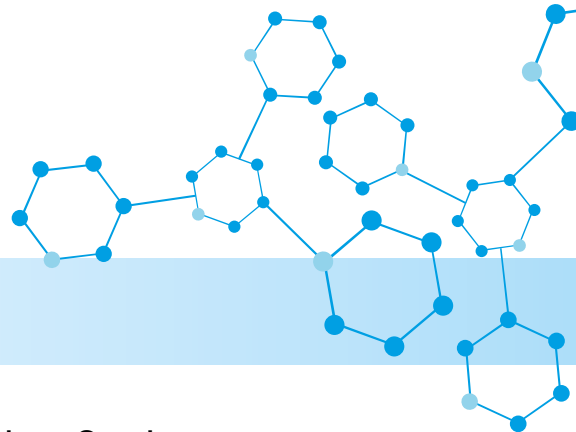

## SPECIFICATIONS

Construction Case	304 Stainless Steel
Construction Tank	304 Stainless Steel
Control	Digital LCD
Overall Dimensions	550 x 490 x 365mm (WxDxH)
Tank Dimensions	505 x 300 x 200mm (WxDxH)
Internal Basket Dimensions	465 x 265 x 160mm (WxDxH)
Maximum Capacity	28Litres
Working Capacity	25Litres
Ultrasonic Power	400W
Ultrasonic Power (per litre)	16W/Litre
Operating Frequency	32 – 38kHz <i>Frequency Leap</i>
Maximum Heating Capacity	Ambient - 70°C
Heater Power	650W
Timer	0 – 99 minutes
Mains Connection	230VAC
Drain Outlet	1/2" BSP valve*
Lid	ABS Plastic
Basket	304 Stainless Steel
Gross Weight (inc. lid, basket & packaging)	18.5kg
Net Weight	17kg
Warranty	3 Years

\*Drain is located at bottom right at the back of the machine

For more information on these products contact: